


## \* MULTICONNECTOR COLD WELDED

⚠ use special WOODFIELD ETO-tape in all cases/ see photo  
özel WOODFIELD Gasband kullanınız / resme bakınız

Gasbandın yarısı PVC üzerine diğer yarısı bezin üzerine gelmeli.  
Glue Gasband 1/2 on Towel and tube

Gasbandı sağdan sola ilerleyerek kademe kademe sarın  
coil Gasband from left to right side around the tube

⚠ Tyvekband see photo  
Tyvekband için resme bakınız!

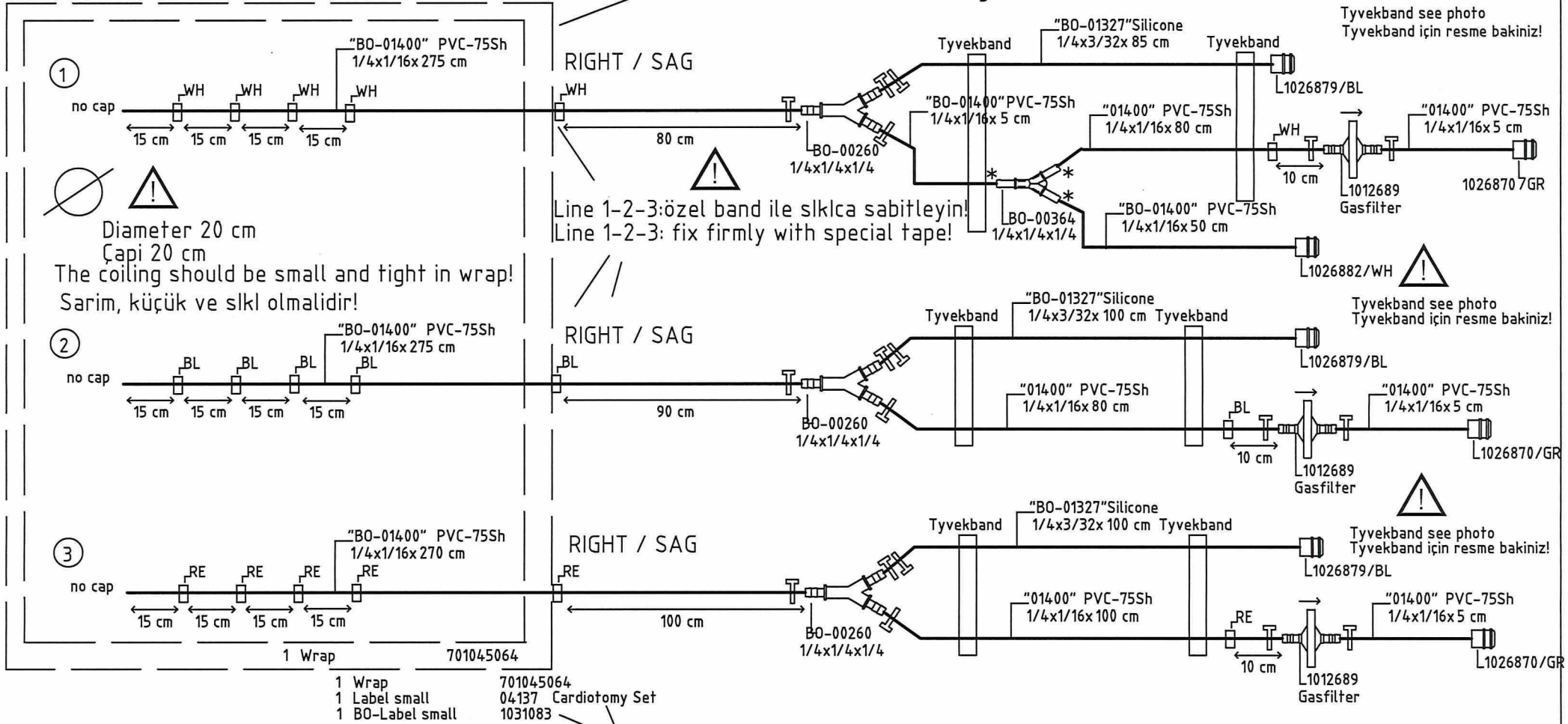

Packaging: ⚠ 04137 ve 701031083 BO-label small etiketini 701045064 wrap üzerine gasband ile sabitleyiniz. <sup>Gastape</sup>  
fix label 04137 and 701031083 BO-label small on 701045064 wrap with gasband.



notice folding - technique /Paketleme teknigine dikkat ediniz!  
notice postion Gasband / Gasbandin pozisyonuna dikkat ediniz!  
SEE PHOTO! / RESME BAKINIZ!

Country NLD City Rotterdam  
Institution Erasmus Medisch Centrum

Circuit Type Adult Pack OPEN  
with HMOD 71000

**MAQUET**  
GETINGE GROUP

TUBING SET


This document contains confidential and proprietary information  
All design, reproduction, manufacturing, use and sale rights in such information are reserved

|                                                |                  |                             |
|------------------------------------------------|------------------|-----------------------------|
| Article-No.<br>701041635                       | Rev.-No.<br>33   | Catalog-No.<br>BO-HQV 51601 |
| FB-713 / Revision 00<br>Released on 2021-01-28 | Date<br>21-12-06 | Drawn by<br>N.Tuncer        |
|                                                |                  | Page<br>1 / 5               |

# \* MULTICONNECTOR COLD WELDED

! use special WOODFIELD ETO-tape in all cases/ see photo  
! özel WOODFIELD Gasband kullanınız / resme bakınız

! Gasbandın yarısı PVC üzerine diğer yarısı bezin üzerine gelmeli.  
Glue Gasband 1/2 on Towel and tube  
Gasbandı sağdan sola ilerleyerek kademe kademe sarın  
coil Gasband from left to right side around the tube

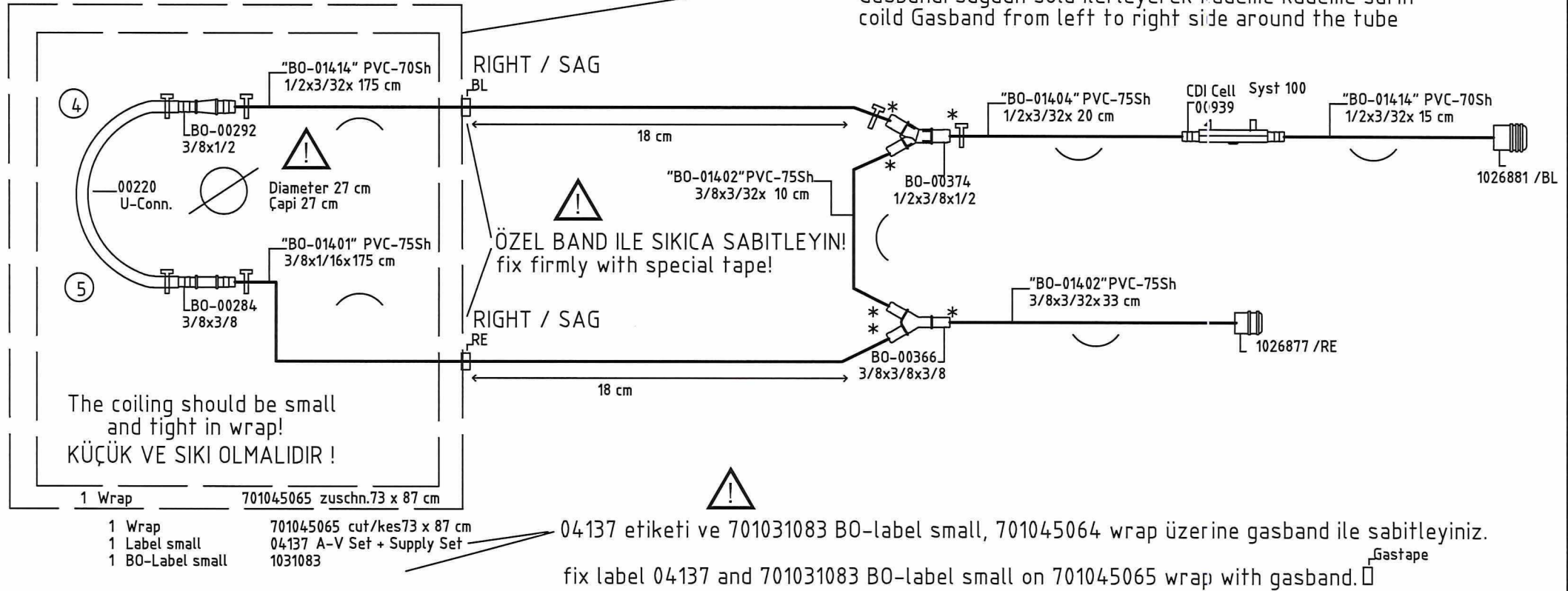

## Packaging:

! notice folding - technique /Paketleme teknigine dikkat ediniz!  
! notice postion Gasband / Gasbandin pozisyonuna dikkat ediniz!  
SEE PHOTO! / RESME BAKINIZ!

Country NLD City Rotterdam  
Institution Erasmus Medisch Centrum

Circuit Type  
Adult Pack OPEN  
with HMOD 71000

# MAQUET

## GETINGE GROUP

### TUBING SET




This document contains confidential and proprietary information  
All design, reproduction, manufacturing, use and sale rights in such information are reserved

|                          |                 |                             |
|--------------------------|-----------------|-----------------------------|
| Article-No.<br>701041635 | Rev.-No.:<br>33 | Catalog-No.<br>BO-HQV 51601 |
|--------------------------|-----------------|-----------------------------|

|                                               |                  |                      |               |
|-----------------------------------------------|------------------|----------------------|---------------|
| FB-71 / Revision 00<br>Released on 2021-01-28 | Date<br>21-12-06 | Drawn by<br>N.Tuncer | Page<br>2 / 5 |
|-----------------------------------------------|------------------|----------------------|---------------|

**\* MULTICONNECTOR COLD WELDED**

TOP / ÜST

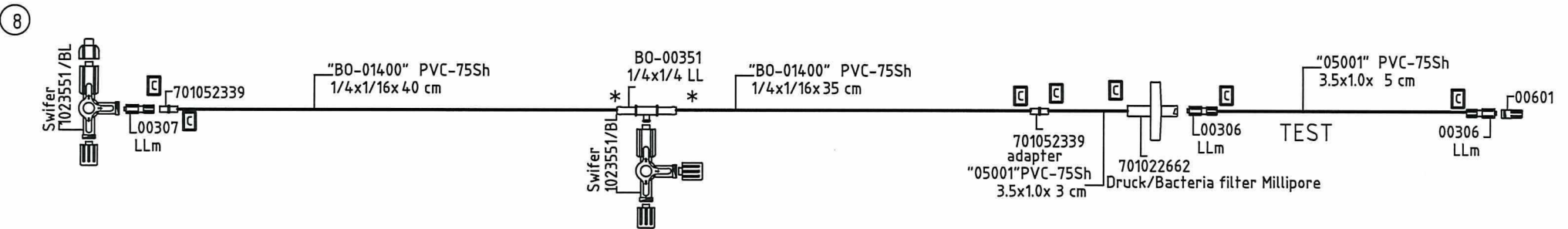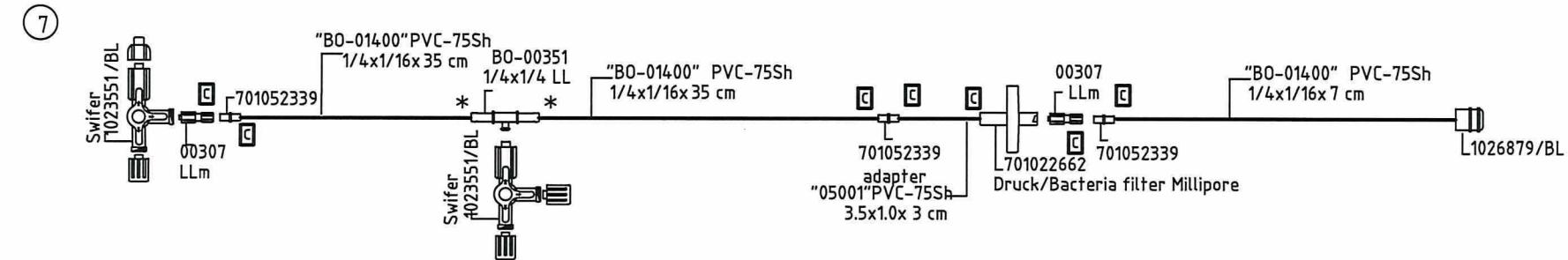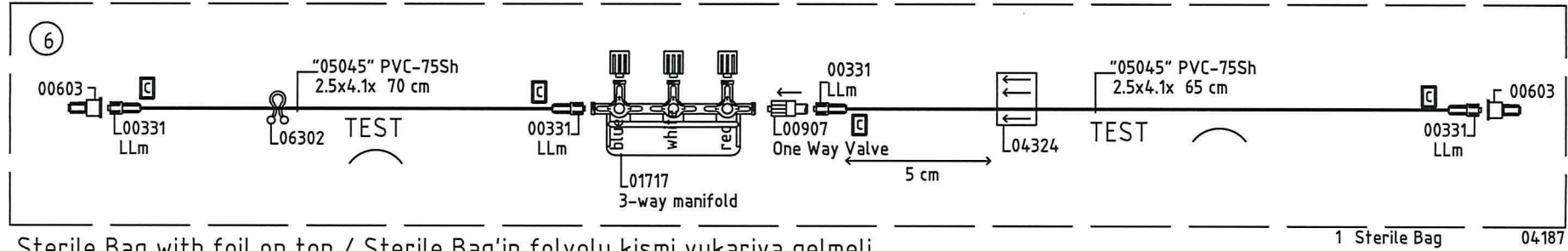

Packaging:

Country  
NLD  
Institution  
Erasmus Medisch Centrum

City  
Rotterdam

Circuit Type  
Adult Pack OPEN  
with HMOD 71000

**MAQUET**  
GETINGE GROUP  
TUBING SET

This document contains confidential and proprietary information  
All design, reproduction, manufacturing, use and sale rights in such information are reserved

|                          |                 |                             |
|--------------------------|-----------------|-----------------------------|
| Article-No.<br>701041635 | Rev.-No.:<br>33 | Catalog-No.<br>BO-HQV 51601 |
|--------------------------|-----------------|-----------------------------|

|                                                |                  |                      |               |
|------------------------------------------------|------------------|----------------------|---------------|
| FB-713 / Revision 00<br>Released on 2021-01-28 | Date<br>21-12-06 | Drawn by<br>N.Tuncer | Page<br>3 / 5 |
|------------------------------------------------|------------------|----------------------|---------------|

⑨




⚠ 701031083 BO-label small, 701045064 wrap üzerine gasband ile sabitleyiniz.  
fix BO-label small 701031083 on 701045065 wrap with gasband. <sup>Gas tape</sup> □

 packaging as Cardiotomy set loop!!!  
Cardiotomy set gibi paketleyiniz!!!

Circuit Type  
Adult Pack OPEN  
with HMOD 71000

|                                                  |                  |                      |               |
|--------------------------------------------------|------------------|----------------------|---------------|
| FB-71(i) / Revision 00<br>Released on 2021-01-28 | Date<br>21-12-06 | Drawn by<br>N.Tuncer | Page<br>4 / 5 |
|--------------------------------------------------|------------------|----------------------|---------------|

# \* MULTICONNECTOR COLD WELDED

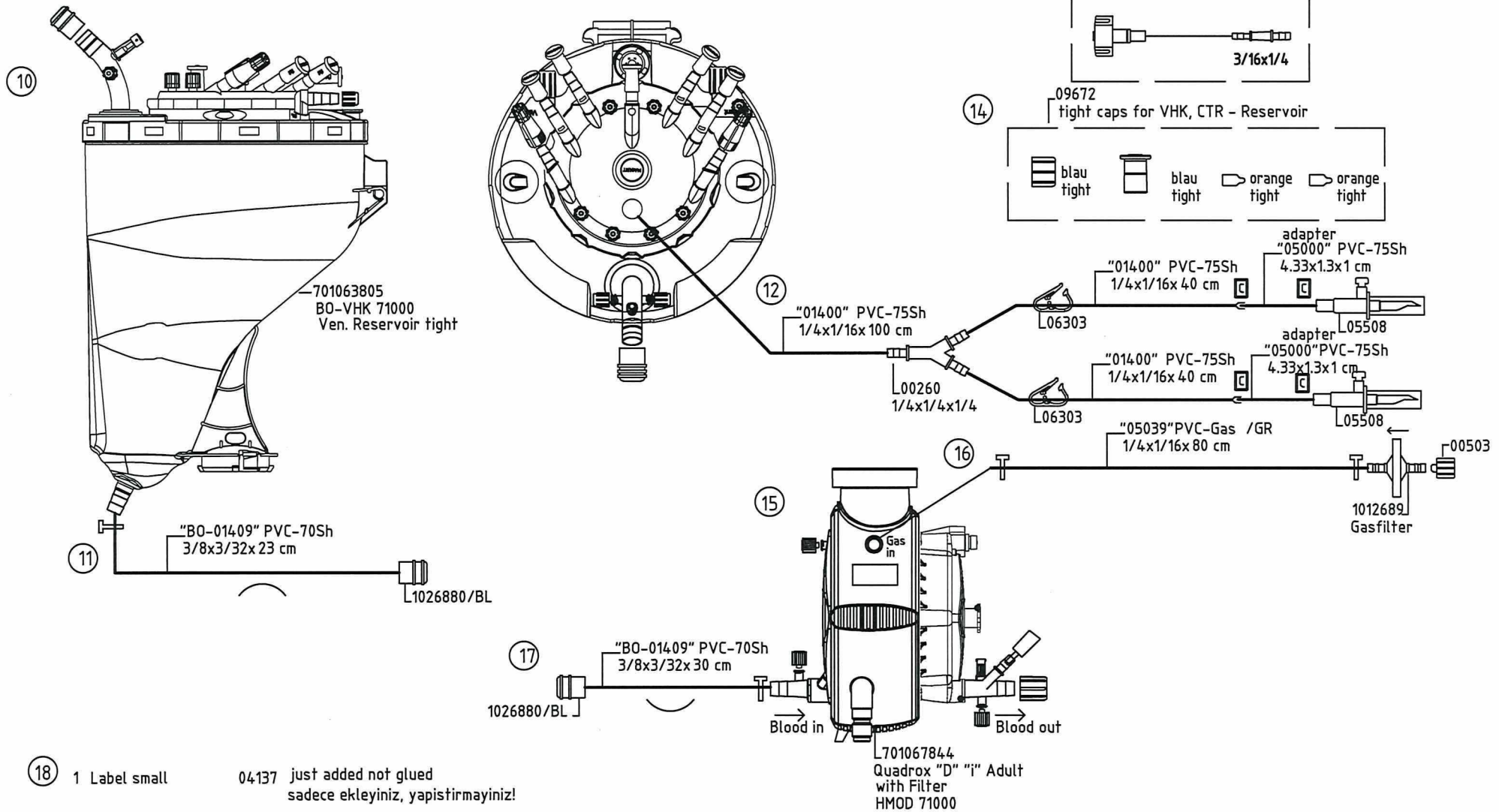

## Packaging:

|                     |           |                  |
|---------------------|-----------|------------------|
| 1 Tray (80x60x26cm) | 07535     | inserti kesiniz. |
| 1 Insert VHK 71000  | 701063333 | cut insert       |
| 1 Velstrap long     | 1032804   |                  |
| 1 Inlay Quadrox-i   | 701063334 |                  |
| 1 Velstrap short    | 1032803   |                  |
| 1 Tyvekcover        | 07560     |                  |
| 2 Label big         | 04148     |                  |
| 2 BO-Label middle   | 1031082   |                  |
| 1 Quality Pass      | 1028298   |                  |
| 1 Lot Label, grey   | 1020612   |                  |

|                       |           |
|-----------------------|-----------|
| 1 Label big JCTS      | 1030276   |
| 1 Paper between       | 07528     |
| 1 Carton              | 07527     |
| 1 IFU HLM Pack        | 04248     |
| 1 IFU VHK 70000/71000 | 701062629 |
| 1 IFU HMOD 71000      | 701067798 |
| 1 IFU Softline        | 1031075   |

Country  
NLD

City  
Rotterdam

Institution

Erasmus Medisch Centrum

Circuit Type  
Adult Pack OPEN  
with HMOD 71000

**MAQUET**  
GETINGE GROUP

## TUBING SET

This document contains confidential and proprietary information  
All design, reproduction, manufacturing, use and sale rights in such information are reserved

|                                                |                  |                             |
|------------------------------------------------|------------------|-----------------------------|
| Article-No.<br>701041635                       | Rev.-No.<br>33   | Catalog-No.<br>BO-HQV 51601 |
| FB-713 / Revision 00<br>Released on 2021-01-28 | Date<br>21-12-06 | Drawn by<br>N.Tuncer        |
|                                                |                  | Page<br>5 / 5               |