


*** MULTICONNECTOR COLD WELDED**

⚠ use special WOODFIELD ETO-tape in all cases/ see photo
özel WOODFIELD Gasband kullanınız / resme bakınız



Gasbandin yarisi PVC üzerine diger yarisi bezin üzerine gelmeli.
Glue Gasband 1/2 on Towel and tube

Gasbandi sagdan sola ilerleyerek kademe kademe sarin
coil Gasband from left to right side around the tube

⚠ Tyvekband see photo
Tyvekband için resme bakınız!

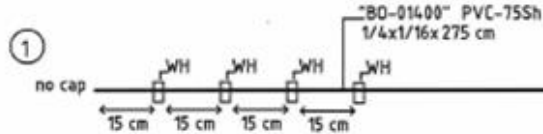

⚠ Diameter 20 cm
Çapı 20 cm
The coiling should be small and tight in wrap!
Sarım, küçük ve sıkı olmalıdır!

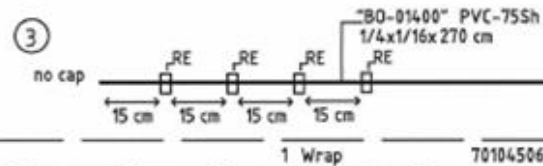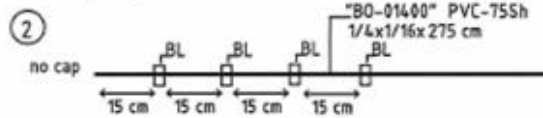


1 Wrap 701045064
1 Label small 04137
1 BO-Label small 1031083

RIGHT / SAG

Line 1-2-3: özel band ile sıkıca sabitleyin!
Line 1-2-3: fix firmly with special tape!

RIGHT / SAG

RIGHT / SAG

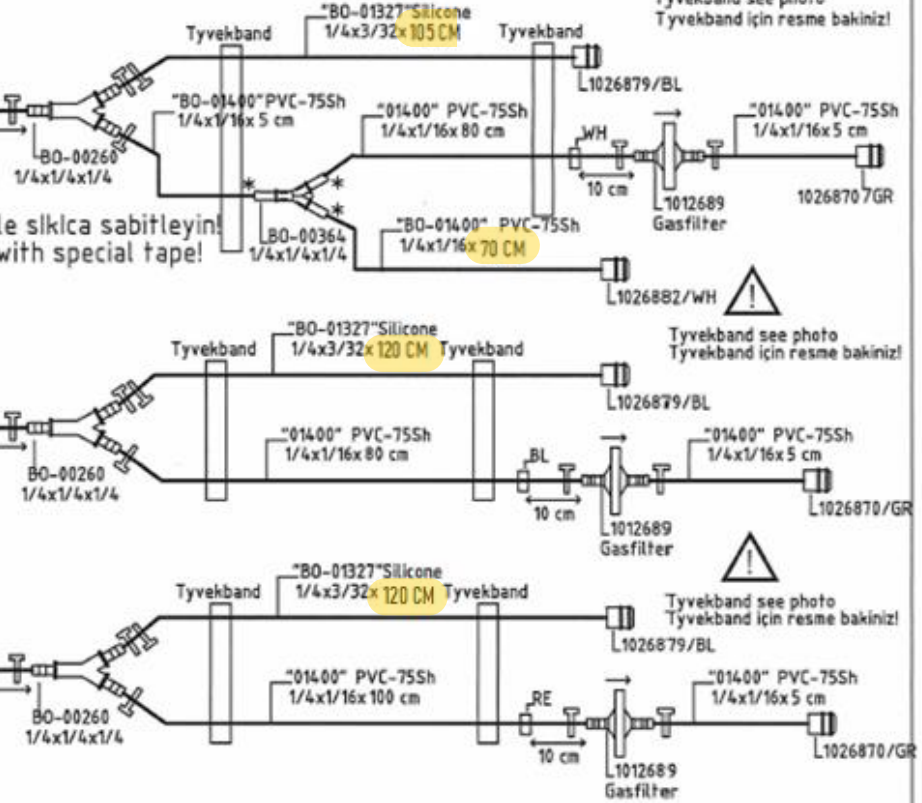

⚠ Tyvekband see photo
Tyvekband için resme bakınız!

⚠ Tyvekband see photo
Tyvekband için resme bakınız!

Packaging: ⚠ 04137 ve 701031083 BO-label small etiketini, 701045064 wrap
üzerine gasband ile sabitleyiniz. Gasband
fix label 04137 and 701031083 BO-label small on 701045064 wrap with gasband.

⚠ notice folding - technique / Paketleme teknigine dikkat ediniz!
notice postion Gasband / Gasbandin pozisyonuna dikkat ediniz!
SEE PHOTO! / RESME BAKINIZ!

Country NLD City Rotterdam
Institution Erasmus Medisch Centrum

Circuit Type
Adult Pack OPEN
with HMOD 71000

MAQUET
GETINGE GROUP

TUBING SET


This document contains confidential and proprietary information.
All design, reproduction, manufacturing, use and sale rights in each information are reserved.

Article-No. 701041635	Rev.-No. 33	Catalog-No. BO-HQV 51601
FB-713 / Revision 00 Released on 2025-05-20	Date 21-12-06	Drawn by N.Tuncer
		Page 1 / 5

* MULTICONNECTOR COLD WELDED

! use special WOODFIELD ETO-tape in all cases/ see photo
! özel WOODFIELD Gasband kullanınız / resme bakınız

! Gasbandın yarısı PVC üzerine diğer yarısı bezin üzerine gelmeli.
Glue Gasband 1/2 on Towel and tube
Gasbandı sağdan sola ilerleyerek kademe kademe sarın
coil Gasband from left to right side around the tube

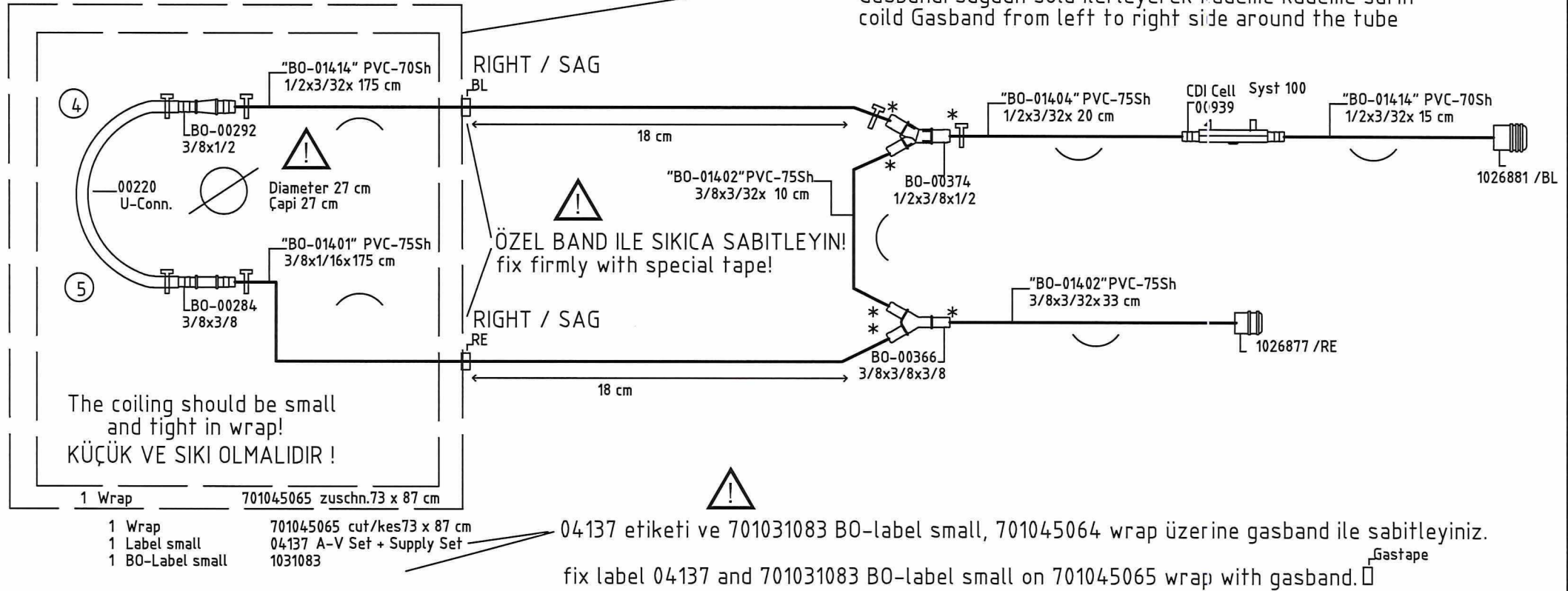

Packaging:

! notice folding - technique /Paketleme teknigine dikkat ediniz!
! notice postion Gasband / Gasbandin pozisyonuna dikkat ediniz!
SEE PHOTO! / RESME BAKINIZ!

Country NLD City Rotterdam
Institution Erasmus Medisch Centrum

Circuit Type Adult Pack OPEN
with HMOD 71000

MAQUET

GETINGE GROUP

TUBING SET




This document contains confidential and proprietary information
All design, reproduction, manufacturing, use and sale rights in such information are reserved

Article-No. 701041635 Rev.-No.: 33 Catalog-No. BO-HQV 51601

FB-71 / Revision 00 Released on 2021-01-28 Date 21-12-06 Drawn by N.Tuncer Page 2 / 5

*** MULTICONNECTOR COLD WELDED**

TOP / ÜST

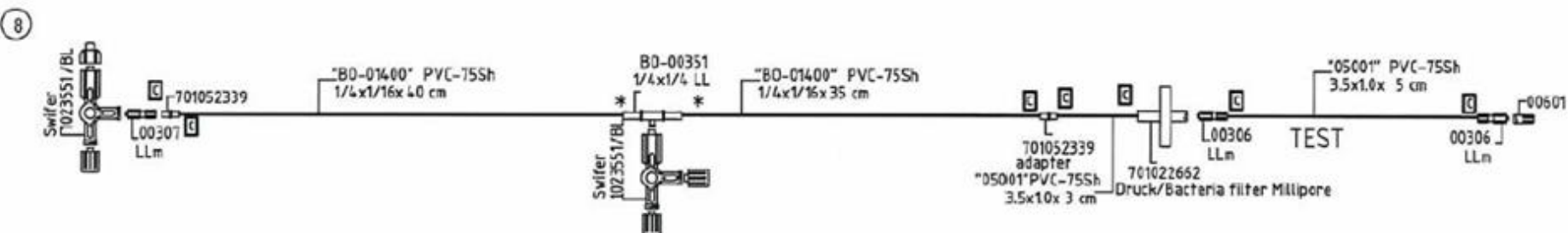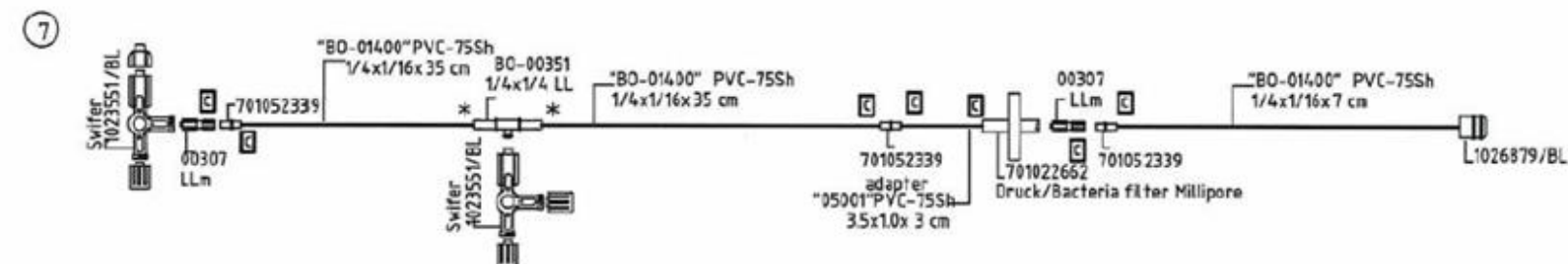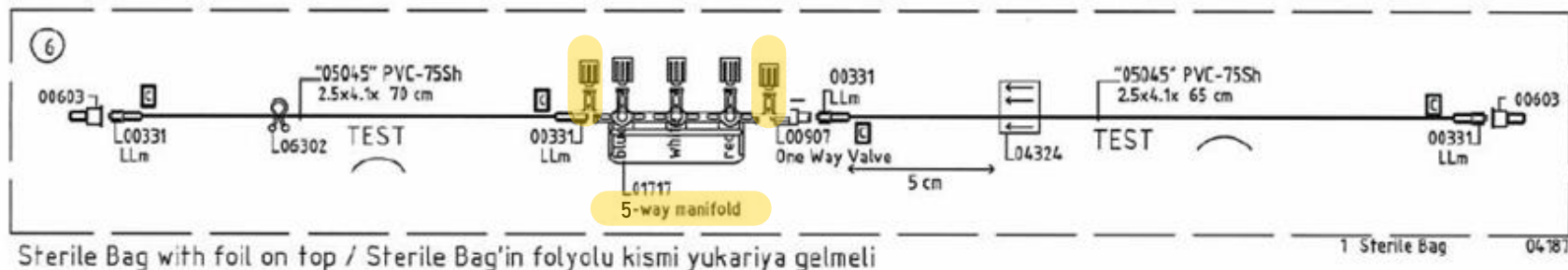

Packaging:

Country
NLD
Institution
Erasmus Medisch Centrum

Circuit Type
Adult Pack OPEN
with HMOD 71000

MAQUET
GETINGE GROUP

TUBING SET

This document contains confidential and proprietary information
All design, reproduction, manufacturing use and sale rights in such information are reserved

Article-No. 70104.1635	Rev.-No. 33	Catalog-No. BO-HQV 51601
---------------------------	----------------	-----------------------------

FB-713 / Revision 10 Released on 2021-01-28	Date 21-12-06	Drawn by N.Tuncer	Page 3 / 5
--	------------------	----------------------	---------------

⑨




⚠ 701031083 BO-label small, 701045064 wrap üzerine gasband ile sabitleyiniz.
fix BO-label small 701031083 on 701045065 wrap with gasband. ^{Gas tape} □

 packaging as Cardiotomy set loop!!!
Cardiotomy set gibi paketleyiniz!!!

Circuit Type
Adult Pack OPEN
with HMOD 71000

FB-711 / Revision 00 Released on 2021-01-28	Date 21-12-06	Drawn by N.Tuncer	Page 4 / 5
--	------------------	----------------------	---------------

*** MULTICONNECTOR COLD WELDED**

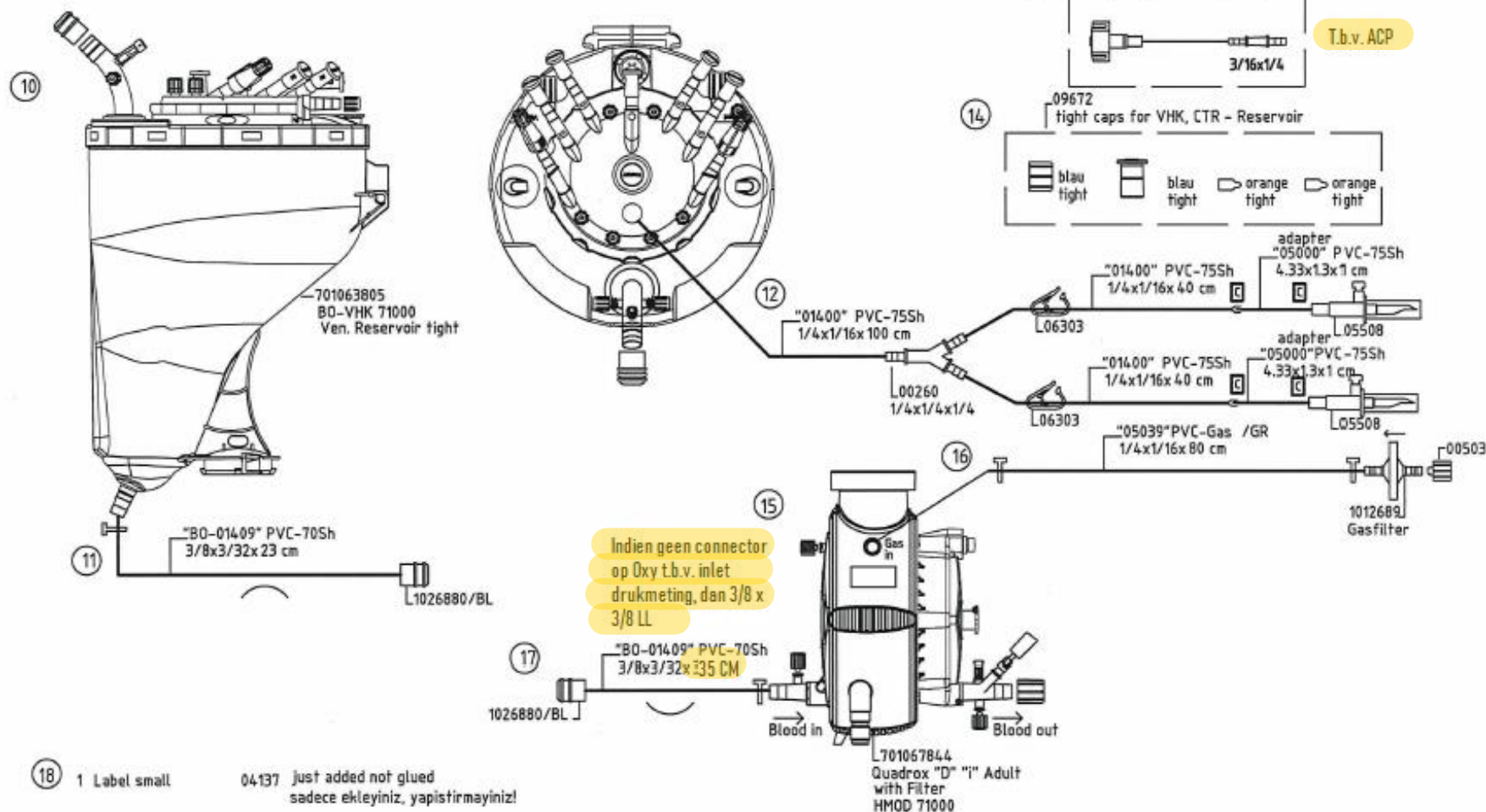

Packaging:

1 Tray (80x60x26cm)	07535	
1 Insert VHK 71000	701063333	inserti kesiniz. cut insert
1 Velstrap long	1032804	
1 Inlay Quadrox-i	701063334	
1 Velstrap short	1032803	
1 Tyvekcover	07560	
2 Label big	04148	
2 BO-Label middle	1031082	
1 Quality Pass	1028298	
1 Lot Label, grey	1020612	

1 Label big JCTS	1030276
1 Paper between	07528
1 Carton	07527
1 IFU HLM Pack	04248
1 IFU VHK 70000/71000	701062629
1 IFU HMOD 71000	701067798
1 IFU Softline	1031075

Country	City
NLD	Rotterdam
Institution	
Erasmus Medisch Centrum	
Circuit Type	
Adult Pack OPEN with HMOD 71000	

MAQUET
GETINGE GROUP

TUBING SET

This document contains confidential and proprietary information.
All design, reproduction, manufacturing, use and sale rights in such information are reserved.

Article-No.	Rev.-No.	Catalog-No.
701041635	33	BO-HQV 51601
FB-T13 / Revision 09 Released on 2021-01-28	Date 21-12-06	Drawn by N.Tuncer
		Page 5 / 5